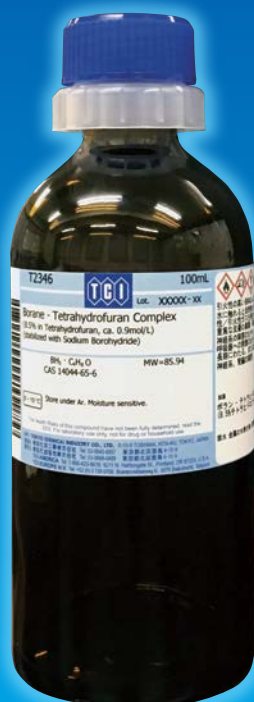


DualSeal

Airtight-Double-Cap Bottled Reagents



Volatile Solutions

Metallic Salt Solutions

Dehydrated Solvents

Organometallic Reagents

Other Reagents

DualSeal

Airtight-Double-Cap Bottled Reagents

TCI introduced our newly developed double cap system "DualSeal" for moisture/oxygen-sensitive products, allowing you to keep them in good condition until the last drop. We will continuously increase our usage of DualSeal across our product portfolio.



Features of DualSeal

- Highly airtight double cap structure.
- Air-sensitive reagents can be safely dispensed without exposure to air.
- Even after piercing the septum cap, the PTFE sheet on the blue outer cap provides sealing protection.
- Easy to dispose of caps and bottles separately.



Outer cap (Blue)

DualSeal specification

DualSeal consists of two parts: the blue outer cap and the white septum cap, both of which can be screwed on and off.

No additional sealing is required after piercing the septum cap with a needle. Just screw the outer cap back in place.

The outer cap has a convex structure on the inside. By filling the space where air and moisture stay, the material is protected from moisture and oxygen even after piercing the septum cap with a needle.

Caution

Do not dispose of the inlaid PTFE sheet! Keep it as it is during use! This PTFE sheet acts as a second layer air seal.



Septum cap (white)

The white septum cap has a wide septum surface for ease of use, and features two layers of rubber and a highly chemically resistant PTFE seal. The septum cap body is made of polypropylene and contains a screw thread allowing for easy removal from the bottle. The septum cap is closed with high torque to ensure an airtight seal. Open the septum cap only when all of the liquid has been used up and you want to prepare for the disposal of the bottle.

Caution

Do not place anything on the septum cap! This will significantly reduce seal quality. The cap and bottle can be separated for disposal. It is easy to dispose of caps and bottles separately. Highly reactive reagents may be residing inside the cap. Please take the necessary precautions to avoid accidents due to exposure to oxygen or moisture. Dispose of the bottle and the caps separately after ensuring that no chemical residue is left behind.

How to use DualSeal

In case using a needle (1):

Solvents except for Halogenated hydrocarbon solvents

1. Clamp and secure the reagent bottle before opening.
2. Carefully unscrew the blue outer cap only. Place the cap near the bottle while in use.
3. To prevent air from entering the container, insert a needle with a balloon filled with an inert gas such as argon or nitrogen. Then insert a syringe needle through the septum surface of the septum cap.
4. Fill the syringe with the required amount of liquid.
5. Remove the syringe, inject the liquid into your reaction vessel and safely dispose of the needle.
6. Take the outer cap and screw it tightly back in place.

*Repeated use can lead to increasing the number of holes or increasing the size of existing holes and will over time lead to an increase in air leakage. To prevent deterioration of the septum, reduce the number of injections as much as possible, or purchase a smaller sized bottle.



In case using a needle (2):

Halogenated hydrocarbon solvents such as dichloromethane

1. Clamp and secure the reagent bottle before opening.
2. Carefully unscrew the blue outer cap only. Place the cap near the bottle while in use.
3. Attach a needle to the PTFE tube for liquid delivery connected to the reaction vessel filled with the inert gas. Puncture the septum with the needle so that the needle tip reaches the bottom of the reagent bottle.
4. A syringe filled with an inert gas or a needle connected to an inert gas pump punctures the septum, and the inert gas is sent into the void of the reagent bottle to send the liquid.
5. Remove the needle of the liquid feeding tube first, and then remove the needle that was feeding the inert gas.
6. Take the outer cap and screw it tightly back in place. However, it is recommended that the minimum number of removals be made, as the holes will be degraded by the vapor of the halogenated hydrocarbon solvent.

How to remove the septum cap

Open the septum cap only when all of the liquid has been used up and you want to prepare for the disposal of the bottle.

1. Clamp and secure the bottle before opening.
2. Open the septum cap by unscrewing. (The septum cap is tightly closed with high torque to ensure airtightness. Using tools such as water pump pliers is recommended.)

*Do not remove the outer cap when you remove the septum cap. Take extra care to avoid any spillage of inner liquid.

Notes on bottle disposal

The cap and bottle can be separated for disposal.

- Remove DualSeal by referring to "How to remove the septum cap".
- Highly reactive reagents may be residing inside the cap. Please take the necessary precautions to avoid accidents due to exposure to oxygen or moisture.
- Dispose of the bottle and the caps separately after ensuring that no chemical residue is left behind.

DualSeal Sealability Test: Moisture Analysis

In order to evaluate the sealability of DualSeal, we periodically measure and monitor the moisture increment by Karl Fischer method after piercing the septum of a 500 mL bottle by a needle.

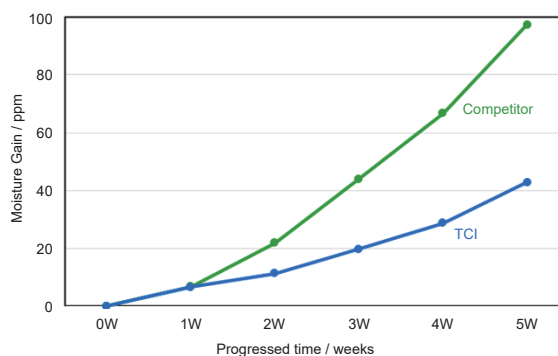
Measurement condition

500 mL anhydrous tetrahydrofuran (TCI product number: T2394) was compared with a competitor's 500 mL anhydrous tetrahydrofuran.

The septum moiety was pierced using an 18 gauge needle at 4 different positions every week (in total, 20 different positions pierced in 5 weeks).

After piercing, the sample solvent was taken and the water content was measured by the Karl Fischer method.

After sampling, the septum was sealed by an outer cap with 1.5 Nm torque and the sealed bottle was stored in a closed environment at ca. 24 °C, ca. 75 % RH (relative humidity).



Periodical measurement to monitor moisture increment every week (5 weeks in total) showed a remarkable difference in the water content (in ppm) between TCI and the competitor's samples. This result indicates that the PTFE sheet inside the outer cap can tightly seal the septum and the outer cap.

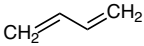
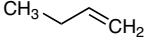
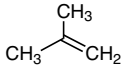
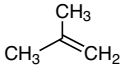
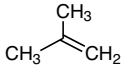
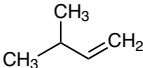
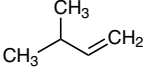
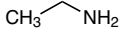
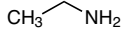
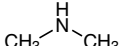
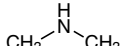
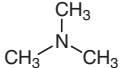
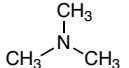
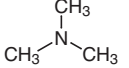
As a result, use of the DualSeal can maintain the quality of product in a sealed bottle for a long time.

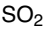

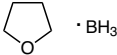






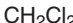
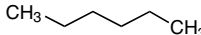

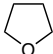


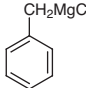


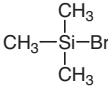
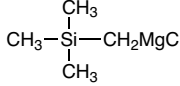
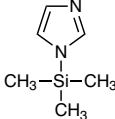
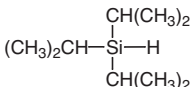
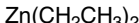
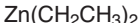
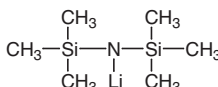
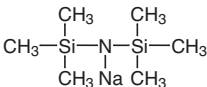
* Not available for smaller than the 100 mL or 100 g size bottles.

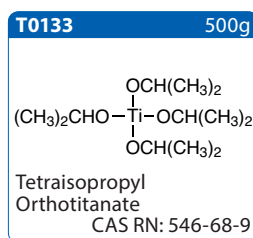
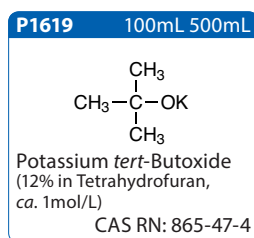
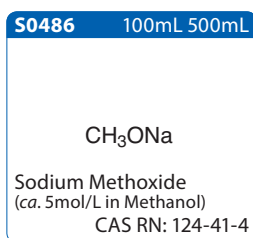
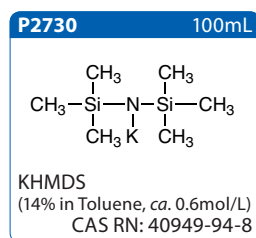
* We are not selling the DualSeal cap itself.

* Since the septum cap part of DualSeal contains butyl rubbers, it is not durable enough for halogenated hydrocarbon solvents such as dichloromethane. Accordingly, the minimum number of needle punctures is recommended, as it is possibly degraded by the vapor of the solvent once a hole is made.

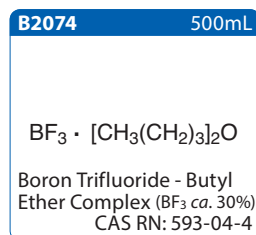
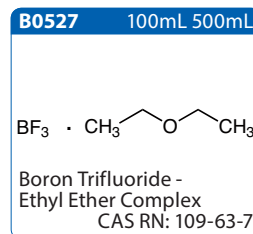
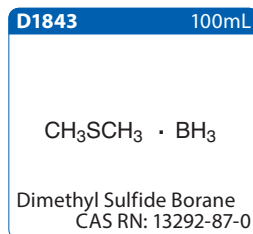
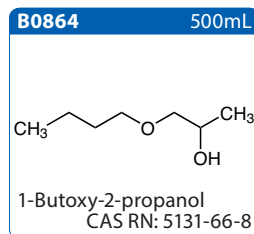
Volatile Solutions

B4358 100mL 500mL				B4359 100mL 500mL		B4410 500mL	
		1,3-Butadiene (ca. 15% in Hexane) CAS RN: 106-99-0		1,3-Butadiene (ca. 15% in Toluene) CAS RN: 106-99-0		1-Butene (ca. 10% in Hexane) CAS RN: 106-98-9	
B4411 100mL		I0909 100mL 500mL		I0910 100mL		I0911 100mL 500mL	
							
1-Butene (ca. 10% in Toluene) CAS RN: 106-98-9		Isobutene (ca. 8% in Dichloromethane) CAS RN: 115-11-7		Isobutene (ca. 10% in Isopropyl Ether) CAS RN: 115-11-7		Isobutene (ca. 15% in Tetrahydrofuran) CAS RN: 115-11-7	
M2563 100mL 500mL		M2565 100mL		P2295 100mL 500mL		P2847 100mL 500mL	
				$\text{CH}_3-\text{C}\equiv\text{CH}$		$\text{CH}_2=\text{C}=\text{CH}_2$	
Isoamylene (ca. 15% in Dichloromethane, ca. 2.5mol/L) CAS RN: 563-45-1		Isoamylene (ca. 12.5% in Tetrahydrofuran, ca. 1.5mol/L) CAS RN: 563-45-1		Propyne (ca. 5% in Tetrahydrofuran, ca. 1mol/L) CAS RN: 74-99-7		Allene (ca. 2% in Tetrahydrofuran, ca. 0.4 mol/L) CAS RN: 463-49-0	
M2813 100mL		T3957 100mL		T3958 100mL		C2882 100mL	
CH_3Cl		CF_3I		CF_3I		$\text{CH}_3\text{CH}_2\text{Cl}$	
Methyl Chloride (ca. 5.7% in Tetrahydrofuran, ca. 1mol/L) CAS RN: 74-87-3		Trifluoroiodomethane (ca. 10% in Tetrahydrofuran, ca. 0.5mol/L) CAS RN: 2314-97-8		Trifluoroiodomethane (ca. 10% in Dimethyl Sulfoxide, ca. 0.6mol/L) CAS RN: 2314-97-8		Chloroethane (ca. 17% in Ethyl Ether, ca. 2.0mol/L) CAS RN: 75-00-3	
V0126 100mL		V0127 100mL		M1016 500mL		M2323 100mL	
$\text{CH}_2=\text{CHBr}$		$\text{CH}_2=\text{CHBr}$		CH_3NH_2		CH_3NH_2	
Vinyl Bromide (ca. 14% in Ethyl Ether, ca. 1.0mol/L) CAS RN: 593-60-2		Vinyl Bromide (ca. 12% in Tetrahydrofuran, ca. 1.0mol/L) CAS RN: 593-60-2		Methylamine (40% in Methanol, ca. 9.8mol/L) CAS RN: 74-89-5		Methylamine (ca. 9% in Ethanol, ca. 2mol/L) CAS RN: 74-89-5	
E0531 100mL		E0817 100mL		D3948 100mL		D4198 100mL	
							
Ethylamine (30-40% in Methanol) CAS RN: 75-04-7		Ethylamine (30-40% in Ethanol) CAS RN: 75-04-7		Dimethylamine (ca. 10% in Tetrahydrofuran, ca. 2mol/L) CAS RN: 124-40-3		Dimethylamine (ca. 11% in Alcohol, ca. 2mol/L) CAS RN: 124-40-3	
T2704 100mL		T3567 100mL 500mL		T2268 100mL		A1884 100mL	
						NH_3	
Trimethylamine (ca. 13% in Tetrahydrofuran, ca. 2mol/L) CAS RN: 75-50-3		Trimethylamine (ca. 13% in Acetonitrile, ca. 2mol/L) CAS RN: 75-50-3		Trimethylamine (ca. 25% in Methanol, ca. 3.2mol/L) CAS RN: 75-50-3		Ammonia (ca. 4% in Methanol, ca. 2.0mol/L) CAS RN: 7664-41-7	
A2236 100mL 500mL		A2237 100mL		H1060 500mL		H1061 100mL	
NH_3		NH_3		HCl		HCl	
Ammonia (ca. 4% in Ethanol, ca. 2.0mol/L) CAS RN: 7664-41-7		Ammonia (ca. 4% in Isopropyl Alcohol, ca. 2.0mol/L) CAS RN: 7664-41-7		Hydrogen Chloride (ca. 1mol/L in Ethyl Acetate) CAS RN: 7647-01-0		Hydrogen Chloride (ca. 1mol/L in Ethyl Ether) CAS RN: 7647-01-0	
H1277 500mL		X0041 100mL 500mL		H0959 100mL 500mL			
HCl		Hydrogen Chloride - Methanol Reagent (5-10%) [for Esterification] CAS RN: 7647-01-0		Hydrogen Bromide - Ethanol Reagent (10-20%) [for Esterification] CAS RN: 10035-10-6			
Hydrogen Chloride (ca. 16% in Cyclopentyl Methyl Ether, ca. 4mol/L) CAS RN: 7647-01-0							

U0147 100mL  Sulfur Dioxide (ca. 8% in Tetrahydrofuran, ca. 1.2 mol/L) CAS RN: 7446-09-5	U0148 100mL  Sulfur Dioxide (ca. 2.5% in Dichloromethane, ca. 0.5 mol/L) CAS RN: 7446-09-5	T2346 100mL 500mL  Tetrahydrofuran Borane (8.5% in Tetrahydrofuran, ca. 0.9mol/L) (stabilized with Sodium Borohydride) CAS RN: 14044-65-6			
Metallic Salt Solutions			L0186 100mL  Lithium Borohydride (ca. 4mol/L in Tetrahydrofuran) CAS RN: 16949-15-8	L0222 100mL  Lithium Chloride (2.3% in Tetrahydrofuran, ca. 0.5mol/L) CAS RN: 7447-41-8	T2052 100mL 500mL  Titanium(IV) Chloride (14% in Dichloromethane, ca. 1.0mol/L) CAS RN: 7550-45-0
T2053 100mL  Tin(IV) Chloride (ca. 1.0mol/L in Dichloromethane) CAS RN: 7646-78-8	Z0019 100mL  Zinc Chloride (ca. 7% in Tetrahydrofuran, ca. 0.5mol/L) CAS RN: 7646-85-7	Z0020 100mL  Zinc Chloride (ca. 25% in 2-Methyltetrahydrofuran, ca. 2mol/L) CAS RN: 7646-85-7			
Dehydrated Solvents			D3478 500mL  Dichloromethane Anhydrous (stabilized with 2-Methyl-2-butene) CAS RN: 75-09-2	H1197 500mL  Hexane Anhydrous CAS RN: 110-54-3	D3479 500mL  Diethyl Ether Anhydrous (stabilized with BHT) CAS RN: 60-29-7
T2394 500mL  Tetrahydrofuran Anhydrous (stabilized with BHT) CAS RN: 109-99-9	Organometallic Reagents				
			M0362 250g  Methylmagnesium Bromide (12% in Tetrahydrofuran, ca. 1mol/L) CAS RN: 75-16-1	E0134 250g  Ethylmagnesium Bromide (39% in Ethyl Ether, ca. 3mol/L) CAS RN: 925-90-6	B1933 250g  Benzylmagnesium Chloride (ca. 16% in Tetrahydrofuran, ca. 1mol/L) CAS RN: 6921-34-2
O0240 250g  n-Octylmagnesium Bromide (ca. 22% in Tetrahydrofuran, ca. 1mol/L) CAS RN: 17049-49-9	S0467 100g 500g  Sodium Bis(2-methoxyethoxy)aluminum Dihydride (70% in Toluene, ca. 3.6mol/L) CAS RN: 22722-98-1	B1087 250mL  Bromotrimethylsilane CAS RN: 2857-97-8	T1451 100mL  Trimethylsilylmethylmagnesium Chloride (20% in Ethyl Ether, ca. 1mol/L) CAS RN: 13170-43-9	T0585 100g  N-Trimethylsilylimidazole CAS RN: 18156-74-6	
T1533 100mL  Triisopropylsilane CAS RN: 6485-79-6	D3214 100mL 500mL  Diethylzinc (ca. 15% in Hexane, ca. 1mol/L) CAS RN: 557-20-0	D3902 100mL  Diethylzinc (ca. 15% in Toluene, ca. 1mol/L) CAS RN: 557-20-0	H0915 100mL 500mL  LiHMDS (ca. 26% in Tetrahydrofuran, ca. 1.3mol/L) CAS RN: 4039-32-1	H0894 100mL 500mL  NaHMDS (contains 2-Methyl-2-butene) (38% in Tetrahydrofuran, ca. 1.9mol/L) CAS RN: 1070-89-9	



Other Reagents



* Not available for smaller than the 100 mL or 100 g size bottles.

* We are not selling the DualSeal cap itself.

* Since the septum cap part of DualSeal contains butyl rubbers, it is not durable enough for halogenated hydrocarbon solvents such as dichloromethane. Accordingly, the minimum number of needle punctures is recommended, as it is possibly degraded by the vapor of the solvent once a hole is made.

Ordering and Customer Service

TCI AMERICA

Tel : 800-423-8616 / 503-283-1681
Fax : 888-520-1075 / 503-283-1987
E-mail : Sales-US@TCIchemicals.com

TCI EUROPE N.V.

Tel : +32 (0)3 735 07 00
Fax : +32 (0)3 735 07 01
E-mail : Sales-EU@TCIchemicals.com

TCI Deutschland GmbH

Tel : +49 (0)6196 64053-00
Fax : +49 (0)6196 64053-01
E-mail : Sales-DE@TCIchemicals.com

Tokyo Chemical Industry UK Ltd.

Tel : +44 (0)1865 784560
E-mail : Sales-UK@TCIchemicals.com

TCI Chemicals (India) Pvt. Ltd.

Tel : 1800 425 7889 / 044-2262 0909
Fax : 044-2262 8902
E-mail : Sales-IN@TCIchemicals.com

梯希爱(上海)化成工业发展有限公司

Tel : 800-988-0390 / 021-67121386
Fax : 021-6712-1385
E-mail : Sales-CN@TCIchemicals.com

TOKYO CHEMICAL INDUSTRY CO., LTD.

Tel : +81 (0)3-5640-8878
E-mail : globalbusiness@TCIchemicals.com

Availability, price or specification of the listed products are subject to change without prior notice. Reproduction forbidden without the prior written consent of Tokyo Chemical Industry Co., Ltd.